Model	Go Go Elite Traveller® 2 Platinum (4-Wheel)
Part Number	SC442EMS
Weight Capacity	300 lbs.
Maximum Speed 2	Up to 4.23 mph
Ground Clearance 4	1.34" at motor 3.15" at mid-deck
Turning Radius 4	37"
Overall Length 4	39.37"
Overall Width4	19. 4 5"
Seat-to-Deck Height	16.38 - 17.83"
Seat Height	16.73"
Tires	Front: 2"x 7" Rear: 2.5"x 8"
Range Per Charge ^{2,3}	12 Ah: Up to 6.6 miles at 300 lbs. Up to 8.8 miles at 220 lbs. 18 Ah: Up to 10.4 miles at 300 lbs. Up to 13.8 miles at 220 lbs.
Product Weight	87 lbs. without batteries 106.5 lbs. with 12 Ah batteries 116.6 lbs. with 18 Ah batteries
Heaviest Piece when Disassembled4	39 lbs. (front section)
Standard Seating	Type: Mesh Back Material: Black Vinyl Dimensions: 17"x17" Weight: 23 lbs.
Motor Output	Rear-wheel drive, sealed transaxle, 24-volt DC motor
Brakes	Regenerative and electromechanical
Battery Requirements5	(Two) 12-volt, sealed lead-acid, non-spillable deep-cycle
Battery Weight	12 Ah: 19.5 lbs. (standard) 18 Ah: 29.6 lbs. 2 amp off-board charger
Battery Charger	Lifetime limited on frame;
Warranty	2-year limited on Electronics; 2-year limited on Drivetrain

HCPCS Code: K0899

- 1 Pride® FDA Class II Medical Devices are designed to aid individuals with mobility impairments.
- 2 Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tire condition. These specifications can be subject to a variance of (+/- 10%). 3 Tested in accordance with ANSI/RESNA, WC Vol2, section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.
- 4 Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or 3%).
- 5 Battery weight may vary based on manufacturer.




704-821-7777 www.ShopMobility.com



HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement.

NOTE: All specifications subject to change without notice. The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of \pm 0-3%. Special order seats or non-standard dimensions will not be fold-flat and depth-adjustable.

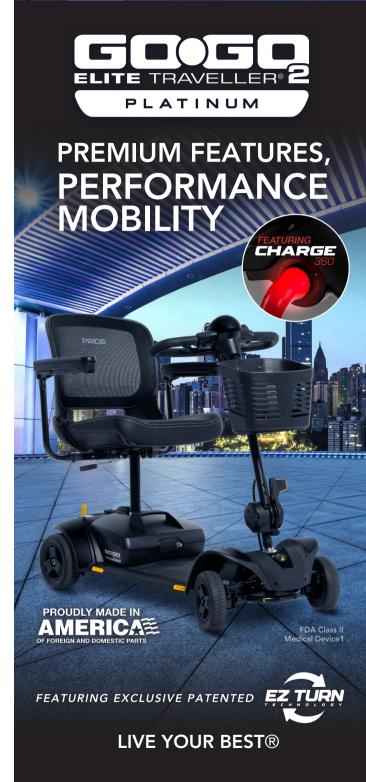



Charge360 technology

FEATURE-RICH, TRAVEL-READY

The Go Go Elite Traveller® 2 Platinum is full of premium features to enhance your mobility experience. With 3-wheel scooter turning ability, magnetic charging, a convenient tiller mechanism, and feather-touch disassembly, the Platinum edition exudes luxury and performance mobility.

- EZ Turn technology for the tight turning of a 3-wheel scooter and the stability of a 4-wheel (turning radius of 37")
- Charge 360 allows you to magnetically charge your scooter, requiring minimal dexterity and no need to line

up any pins

- Angle-adjustable tiller mechanism can be operated with your hand or foot
- Feather-touch disassembly allows the unit to be easily separated into 5 lightweight pieces
- Front frame-mounted seat post for optimal stability
- Removeable footboard basket for additional storage
- 300-lb. weight capacity
- Maximum speed of up to 4.23 mph
- Front basket (standard)

2 sets of removable color-through shrouds included:

Cardinal Red

Black (matte)

Accessories*

- 2.1A XLR USB charger
- Cup holder
- Forearm crutch holder
- Military patches
- Oxygen holder
- Rear basket
- Rear mount 3" double cane/crutch holder
- Rear view mirror

- Saddle bag
- Safety flag
- Single rear mount 3" cane/crutch holder
- Under seat storage
- Walker holder
- Weather cover
- Wishbone crutch holder
- Additional shroud panels

Accessory Shroud Panel Colors

(additional charge applies)



* Only one rear accessory can be used at a time.